


Measurement Report

Surface Texture Measurement

Auto Correlation of Roughness Dataset (Selection)

ZONA_4_1

Lc = 800.000µm

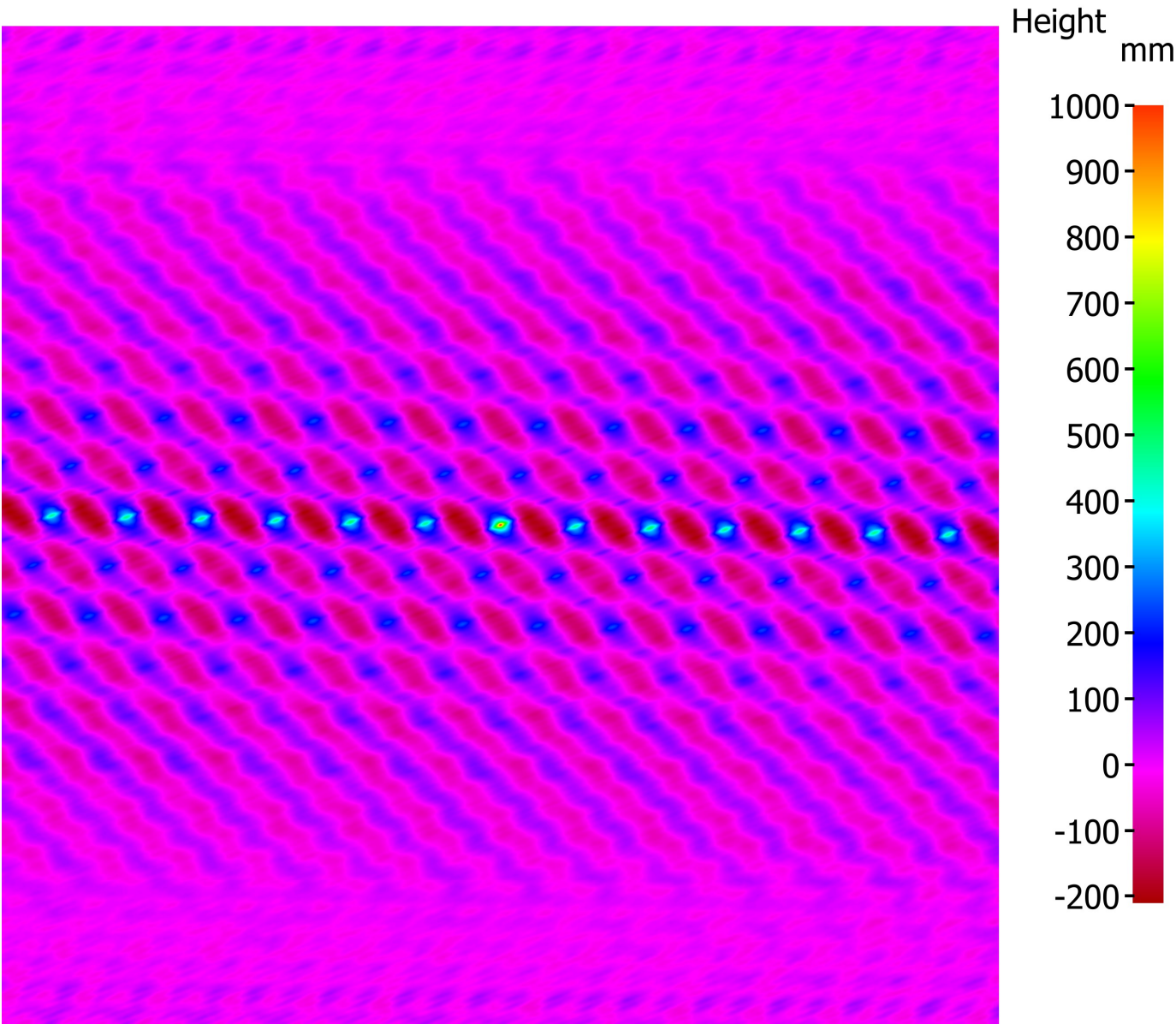

Name	Value	Unit	Description
Sal	40.540	µm	Auto correlation length. Horizontal distance of the auto correlation function (ACF) which has the fastest decay to the value 0.2. Large value: surface dominated by low frequencies. Small value: surface dominated by high frequencies.
Str	0.750		Texture aspect ratio (0 <= Str <= 1). Ratio from the distance with the fastest to the distance with the slowest decay of the ACF to the value 0.2. Str < 0.3: surface has a strong directional structure. Str > 0.5: surface has rather uniform texture.
Std	-1.300	°	Texture direction (°). Derived from the maximum of the angular power spectrum. Std = 90° means a dominant lay parallel to the y-axis.
Stdi	0.615		Texture direction index (0 <= Stdi <= 1). Average value of the angular power spectrum divided by the value of the dominating direction. Smaller values correspond to stronger directional structures.

Name	Value	Unit	Description
angle	n.a.	°	Manually measured angle